# ORIGINAL PAPER

# Milk and serum standard reference materials for monitoring organic contaminants in human samples

Michele M. Schantz · Gauthier Eppe · Jean-François Focant · Coreen Hamilton · N. Alan Heckert · Rebecca M. Heltsley · Dale Hoover · Jennifer M. Keller · Stefan D. Leigh · Donald G. Patterson Jr · Adam L. Pintar · Katherine E. Sharpless · Andreas Sjödin · Wayman E. Turner · Stacy S. Vander Pol · Stephen A. Wise

Received: 10 September 2012 / Revised: 18 October 2012 / Accepted: 23 October 2012 / Published online: 7 November 2012 © Springer-Verlag Berlin Heidelberg (outside the USA) 2012

**Abstract** Four new Standard Reference Materials (SRMs) have been developed to assist in the quality assurance of chemical contaminant measurements required for human biomonitoring studies, SRM 1953 Organic Contaminants in Non-Fortified Human Milk, SRM 1954 Organic Contaminants in Fortified Human Serum, and SRM 1958 Organic Contaminants in Fortified Human Serum, These materials were

developed as part of a collaboration between the National Institute of Standards and Technology (NIST) and the Centers for Disease Control and Prevention (CDC) with both agencies contributing data used in the certification of mass fraction values for a wide range of organic contaminants including polychlorinated biphenyl (PCB) congeners, chlorinated pesticides, polybrominated diphenyl ether (PBDE) congeners, and polychlorinated dibenzo-*p*-dioxin (PCDD) and dibenzo-furan

Contribution of the US Government; not subject to copyright

**Electronic supplementary material** The online version of this article (doi:10.1007/s00216-012-6524-3) contains supplementary material, which is available to authorized users.

M. M. Schantz ( ) · K. E. Sharpless · S. A. Wise Analytical Chemistry Division, National Institute of Standards and Technology (NIST), 100 Bureau Drive, Gaithersburg, MD 20899, USA e-mail: michele.schantz@nist.gov

G. Eppe · J.-F. Focant Mass Spectrometry Laboratory, University of Liège, Allée de la Chimie, 4000 Liège, Belgium

C. Hamilton · D. Hoover AXYS Analytical Services Ltd, 2045 Mills Road West, Sidney, BC V8L 5X2, Canada

N. A. Heckert · S. D. Leigh · A. L. Pintar Statistical Engineering Division, NIST, 100 Bureau Drive, Gaithersburg, MD 20899, USA

R. M. Heltsley · J. M. Keller · S. S. Vander Pol Analytical Chemistry Division, NIST, 331 Fort Johnson Road, Charleston, SC 29412, USA D. G. Patterson Jr · A. Sjödin · W. E. Turner Organic Analytical Toxicology, Division of Laboratory Sciences, National Center for Environmental Health, Centers for Disease Control and Prevention, 4770 Buford Highway, Atlanta, GA 30341, USA

Present Address:
R. M. Heltsley
Aegis Sciences Corporation,
345 Hill Avenue,
Nashville, TN 37210, USA

Present Address:
D. G. Patterson Jr
EnviroSolutions Consulting, Inc.,
4921 Bull Creek Road,
Auburn, GA 30011, USA



(PCDF) congeners. The certified mass fractions of the organic contaminants in unfortified samples, SRM 1953 and SRM 1957, ranged from 12 ng/kg to 2200 ng/kg with the exception of 4,4'-DDE in SRM 1953 at 7400 ng/kg with expanded uncertainties generally <14 %. This agreement suggests that there were no significant biases existing among the multiple methods used for analysis.

**Keywords** Human Serum · Human Milk · SRMs · PCBs · Pesticides PBDEs · PCDFs · PCDDs

## Introduction

Data for contaminants in the human population from biomonitoring studies based on analysis of human serum, milk, adipose tissue, urine, and/or other tissues or fluids [1-8, to cite a few] abound in the literature. These representative studies have been done in the U.S. [1, 2], Canada [3], Belgium [4], Japan [5], Korea [6], Sweden [7], and the Ukraine [8]. The published data are often used by risk assessors to evaluate potential biological effects from the amount of chemicals absorbed into the human body and to identify long-term trends in the population [9]. To compare the data among studies conducted both nationally and internationally, it is important that the accuracy of biomonitoring data be supported by appropriate quality control samples, including the use of reference materials and participation in interlaboratory comparison studies [10]. Risk assessors need to have confidence in the analytical methods that were used to generate the data as well as in the sampling strategy that was used to obtain the samples [10]. As part of the National Health and Nutrition Examination Survey (NHANES) and other biomonitoring programs, the Centers for Disease Control and Prevention (CDC) uses an extensive quality assurance program with as many as 60 % quality control samples included in their analyses [11].

Certified Reference Materials (CRMs), of which Standard Reference Materials (SRMs) are a subset produced by the National Institute of Standards and Technology (NIST), are valuable tools in developing and validating analytical methods to improve quality assurance standards. NIST developed the first serum-matrix SRM for the measurement of organic contaminants in 1985. This material, SRM 1589 Polychlorinated Biphenyls (as Aroclor 1260) in Human Serum [12], was prepared by fortifying human serum with Aroclor 1260 prior to freeze-drying. In 2000 SRM 1589 was replaced with SRM 1589a PCBs, Pesticides and Dioxin/Furans in Human Serum which was prepared using a serum pool collected in 1996 [13] that had been prescreened and found to contain measurable levels of the contaminants of interest. SRM 1589a was not fortified, containing only natural levels of the organic contaminants with mass fractions value assigned for selected polychlorinated biphenyl (PCB) congeners, chlorinated pesticides,

polychlorinated dibenzo-p-dioxin (PCDD) congeners, and polychlorinated dibenzofuran (PCDF) congeners. In 2006 the Certificate of Analysis for SRM 1589a was updated to increase the number of PCB congeners and chlorinated pesticides that were value assigned and also to add mass fraction values for selected polybrominated diphenyl ether (PBDE) congeners [14, 15]. This material found widespread use as a control material for human biomonitoring programs such as the NHANES and is no longer available.

Another matrix of interest in monitoring the exposure of the human population to environmental contaminants is human milk. Human milk is an important exposure route for organic contaminants to the nursing infant [5, 8, 16]. Human milk has approximately 10 times the lipid content compared to human serum; therefore, on a wet-mass basis, serum is not a good surrogate for milk. Some studies, such as Inoue et al. [5] and Hites [16], suggest that similar values of lipophilic contaminants are obtained from serum and breast milk if the measurements are normalized to lipid content. However, the calculation of total lipid varies among laboratories with gravimetric lipid determinations being the traditional method of choice, although some laboratories now use enzymatic analysis of lipid species in serum followed by summation [17]. The differences in lipid content and methods of lipid determinations suggest that quality control samples should match the matrix of interest.

NIST and the CDC Division of Laboratory Sciences have collaborated to develop four new SRMs to meet the expanding needs for the measurement of organic contaminants in human serum and milk. The timing was based not only on the impending depletion of the supply of SRM 1589a but also on the growing need for an SRM representative of more contemporary serum contaminant levels. The goals of the NIST/CDC collaboration were to produce human serum and milk SRMs at current contaminant levels and at elevated levels and to produce a sufficient supply to be available for 10 years. To accomplish these goals, collections of serum (total of 200 L) and milk (total of 100 L) were obtained from blood banks and milk banks, respectively, located around the U.S., combined, and mixed. For each matrix (serum and milk), one SRM is a natural level (SRM 1953 Organic Contaminants in Non-Fortified Human Milk [18] and SRM 1957 Organic Contaminants in Non-Fortified Human Serum [19]), and the other SRM (SRM 1954 Organic Contaminants in Fortified Human Milk [20] and SRM 1958 Organic Contaminants in Fortified Serum [21]) is from the same serum or milk pool to which 169 organic contaminants have been added (See Table S1). The fortified serum and milk SRMs were prepared to represent contaminant concentrations that are approximately ten times higher than the median concentrations found in the U.S. population during the 2003 NHANES study [2].

The preparation and characterization of these four SRMs for PCBs, chlorinated pesticides, PBDEs, PCDDs, and PCDFs are described in this paper. The assignment of mass



fraction values for selected contaminants in these materials involved the combination of results from gas chromatography/mass spectrometry (GC/MS) analyses by NIST and CDC, as well as results from four laboratories, AXYS Analytical Services (Sidney, B.C., Canada), National Institute of Public Health of Québec (Québec, Canada), Stockholm University (Stockholm, Sweden), and University of Liege (Liege, Belgium), which participated in an interlaboratory study using the four SRMs as unknown samples with each laboratory using its normal methods of analysis.

# **Experimental**

Source and preparation of serum and milk SRMs

Plasma used for SRM 1957 and SRM 1958 was acquired in 2004 from various blood banks located around the U.S.: Wilmington and Greenville, NC (10 %); Jacksonville and Orlando, FL (8 %); Jonesboro, AR (12 %); Flagstaff, AZ (15 %); Gallup and Albuquerque, NM (20 %); Memphis, TN (12 %); Portland, ME (16 %); and Carbondale, IL (7 %). Following the precipitation of fibrin and filtration to ensure better homogeneity of the sample, the serum was pooled (approximately 200 L total) and stored at 4 °C. The pool was split into two batches of 100 L each for production of SRM 1957 and SRM 1958. For SRM 1958 approximately 710 mL of a methanol solution containing 169 compounds (Table S1) at varying mass/volume (see below) was added to the 100 L batch of serum, which was then stirred for 4 h. Using a calibrated automatic pipetter, 10.7 mL aliquots of serum were dispensed into 30 mL amber glass vials, and the samples were lyophilized until a stable vacuum and temperature were achieved. Approximately 10,000 vials of each SRM were prepared, and a sales unit of SRM 1957 and SRM 1958 consists of 5 vials of freeze-dried serum. Before use, the serum in each vial must be reconstituted with 10.7 mL of distilled or HPLC-grade water.

The pasteurized milk used for SRM 1953 and SRM 1954 was acquired in 2006 from six milk banks located around the U.S: Florida (4 %), North Carolina (6 %), Iowa (6 %), Delaware (7 %), California (12 %), and Texas (65 %). The SRMs were prepared from milk that had either reached its expiration date or was considered "research milk" that could not be fed to babies. (Research milk includes milk that has too low a caloric content or tested positive for bacteria prior to pasteurization.) The milk was pooled (approximately 100 L total) and was stored at 4 °C. The pool was split into two batches of 50 L each for production of SRM 1953 and SRM 1954. For SRM 1954 approximately 355 mL of a methanol solution containing 169 compounds at varying concentrations (the same solution as added to SRM 1958, see below for more information) was added to the 50 L batch of milk, which was then stirred for 4 h. Using a calibrated automatic pipetter, 5-mL aliquots of milk were dispensed into 10 mL amber glass vials which were then stored at -20 °C. Approximately 5,000 vials of each milk SRM were prepared, and a sales unit of SRM 1953 and SRM 1954 consists of 5 vials of frozen milk.

The methanol spiking solution was prepared by Cambridge Isotope Labs (Andover, MA) and contained 169 compounds at the indicated mass concentration: 38 non and mono ortho-PCBs (50 pg/mL to 500 pg/mL), 4 nonortho-PCBs (0.4 pg/mL to 0.8 pg/mL), 22 chlorinated pesticides including 6 toxaphene congeners (500 pg/mL), 17 PBDEs (500 pg/mL), polybrominated biphenyl (PBB) 153 (500 pg/mL), hexabromocyclododecanes (500 pg/mL), 1,2bis(2,4,6-tribromophenyloxy)ethane (500 pg/mL), hexabromobenzene (500 pg/mL), decabromodiphenylethane (500 pg/mL), 17 PCDDs/PCDFs (0.1 pg/mL to 2.4 pg/mL), 7 polybrominated dibenzodioxins, 10 polybrominated dibenzofurans, and 8 chlorobromo dibenzodioxins and dibenzofurans (0.05 pg/mL), 9 polychlorinated naphthalenes (PCNs) (1 pg/mL), 12 halogenated phenolic compounds, 5 hydroxylated PCBs (500 pg/mL), and 8 other persistent organochlorine compounds (500 pg/mL) including 4 chlorobenzenes, octachlorostyrene, pentachloronitrobenzene, hexachlorocyclopentadiene, and hexachloro-1,3-butadiene. The individual compounds in the spiking solution are listed in Table S1.

Analytical methods The analytical methods used are described in detail in the supplementary information (see Electronic Supplementary Material Figure S1) with a summary here. Multiple methods of analysis were used at NIST for the determination of the PCBs, pesticides, and PBDEs in both the milk and serum, and an additional method was developed for the determination of the hydroxylated compounds in serum only (see Supplemental information). The first method was used for the homogeneity study, analyzing 10 bottles using a liquid-liquid extraction with *n*-hexane: methyl-tert-butyl ether, a sample clean-up using concentrated sulfuric acid followed by silica solid-phase extraction (SPE) fractionation. The final analysis used gas chromatography/mass spectrometry (GC/MS) operated both in the electron ionization (EI) and negative ion chemical ionization (NICI) mode with on-column injections. The GC column used for the EI analysis was a non-polar proprietary phase (Agilent Technologies, Wilmington, DE) while the column used for the NICI analysis was a 50 % (mole fraction) phenyl-substituted methylpolysiloxane phase.

The second method at NIST used focused microwave extraction with 3 mL of 20 % (volume fraction) dichloromethane in hexane as described by Keller et al. [22]. An acidified silica column clean-up followed by an alumina column (5 % deactivated) clean-up was used with final analysis using GC/MS in the EI mode with a 5 % phenyl-substituted



methylpolysiloxane phase. The injections used a programmable temperature vaporization (PTV) inlet. The extracts were also injected on-column using a shorter column (10 m compared to 30 m; note that the short column was used to reduce degradation of the higher brominated PBDEs on-column) connected to the NICI source in the GC/MS.

The methods used by CDC are described in more detail in Patterson and Turner [23] and Sjödin et al. [24, 25]. Gas chromatography/high resolution mass spectrometry (GC/HRMS) with mass resolution of 10,000 and equipped with a 5 % phenyl-substituted methylpolysiloxane phase was used.

The laboratories participating in the interlaboratory study employed the methods commonly used in their laboratories for the determination of these compounds in serum and milk. The methods used by University of Liege are described in more detail in Focant et al. [26].

## Results and discussion

Since the development of the first natural-matrix SRMs for organic contaminants at NIST over 30 years ago, the assignment of certified mass fraction values have been based on the combination of results from two or more independent methods [28]. Currently NIST SRMs for chemical composition are assigned mass fraction values designated certified, reference, or information values based on the number and independence of the analytical methods used, source of the results, and degree of confidence in the accuracy of the assigned value [27, 28].

As described above, results from two to five methods of analysis were used to assign certified and reference mass fraction values for the milk and serum materials. For the PCBs, chlorinated pesticides, and PBDEs, the assigned values were based on results from two or three NIST methods, the CDC method, and the interlaboratory study. The various methods included different extraction procedures, clean up procedures, and GC/MS or GC/HRMS with different stationary phase columns. The PCDD/PCDF congeners were assigned reference values using results from only the CDC method and the interlaboratory study of results from two laboratories. Similarly the reference values for the hydroxylated compounds were assigned based on results from LC/MS/MS (NIST method 3, see Supplemental information) and the interlaboratory study.

The results from each of the five methods are summarized in Electronic Supplementary Material Tables S2 through S5 for SRM 1953, SRM 1954, SRM 1957, and SRM 1958, respectively. The results of the different methods were combined as weighted means [29] to assign the certified and reference mass fraction values. The resulting certified, reference, and information mass fraction values are summarized in

Tables 1, 2, 3, and 4 for the PCB congeners, pesticides, PBDE congeners, and PCDD and PCDF congeners, respectively.

For each of the materials, the homogeneity of the PCBs, pesticides, and PBDEs was assessed at NIST using method 1 described above. An analysis of variance did not show inhomogeneity for a 2.5 g sample of milk or for a 5 g sample of serum [18–21]. Other analytes were treated as though they were homogeneously distributed in the material although homogeneity was not assessed. The long-term stability of the freeze-dried serum materials has been demonstrated for SRM 1589a in which the PCBs and pesticides have been shown to be stable over an 8 year period [15]. Stability assessments are conducted regularly on the freeze-dried serum and frozen milk samples primarily through their use as control samples for other analyses.

For the non-fortified serum and milk, the certified mass fractions of the organic contaminants ranged from 12 ng/kg to 2200 ng/kg (excluding 4,4'-DDE at 7400 ng/kg) with expanded relative uncertainties from approximately 1 % to 14 %, with the exception of PCB 138 in SRM 1957 for which the expanded relative uncertainty is 24 %. Of the 46 certified values in the non-fortified milk and serum SRMs, 26 had relative uncertainties less than 5 % and only 8 had relative uncertainties greater than 10 %. For the fortified serum and milk, the certified mass fractions ranged from 400 ng/kg to 2600 ng/kg (excluding 4,4'-DDE at 8200 ng/ kg) with expanded relative uncertainties from approximately 1 % to 15 %. For the 128 certified values in the fortified milk and serum SRMs, 111 had relative uncertainties less than 10 % with 65 having relative uncertainties less than 5 %. The agreement of the results from the different methods for both the non-fortified and the fortified serum and milk, based on the uncertainties associated with the certified values, are consistent over a wide range in mass fractions (two orders of magnitude) suggesting that no significant biases exist among the methods using different extraction, isolation, and GC/MS columns.

For the natural (non-fortified) serum and milk SRMs, the wet-mass based contaminant levels are higher in the milk compared to the serum due to the higher lipid content of the milk. The milk sample contains a mass fraction of 3.21 %± 0.33 % total extractable lipid as determined by a gravimetric drying step following extraction. Despite using research milk (milk that has too low a caloric content to be fed to infants), the lipid content of the SRM is similar to that found in breast milk from biomonitoring studies [8, 25]. The serum sample contains 0.471 %±0.065 % lipid for SRM 1957 and 0.406 %±0.036 % lipid for SRM 1958 as determined by three methods: a gravimetric drying, the Bligh and Dyer method, and an enzymatic analysis [17]. Comparing the contaminant levels on a lipid-basis, the non-fortified milk sample (SRM 1953) has higher values of the tetra-, penta-, and some hexachlorinated PCB congeners while the



Table 1 Certified, Reference, and *Information* mass fractions (ng/kg milk or reconstituted serum) of selected PCB congeners in SRM 1953 organic contaminants in non-fortified human milk, SRM 1954 organic contaminants in fortified human milk, SRM 1957 organic contaminants in non-fortified human serum, and SRM 1958 organic contaminants in fortified human serum

<sup>a</sup>PCB congeners are numbered according to the scheme proposed by Ballschmiter and Zell [33] and later revised by Schulte and Malisch [34] to conform with IUPAC rules. For the specific congeners mentioned in this table, the Ballschmiter-Zell numbers correspond to those of Schulte and Malisch. When two or more congeners are known to coelute under the GC analysis conditions used, the PCB congener listed first is the major component and the additional congeners may be present as minor components

<sup>b</sup>Values are weighted means of the mass fractions determined by two to five analytical methods [29]. The uncertainty listed with each value is an expanded uncertainty about the mean [29, 35], with coverage factor, k=2 (except k=4 for PCB 180/193 in SRM 1957), calculated by combining a pooled within method variance with a between method variance [29, 36] following the ISO Guide [37, 38]

cValues are the means of results using one analytical technique. The expanded uncertainty, U, is calculated as  $U = ku_c$ , where  $u_c$  is one standard deviation of the analyte mean, and the coverage factor, k, is determined from the Student's t-distribution corresponding to the associated degrees of freedom and 95 % confidence level for each analyte

| PCB congener <sup>a</sup> | SRM 1953                     | SRM 1954                           | SRM 1957                              | SRM 1958                            |
|---------------------------|------------------------------|------------------------------------|---------------------------------------|-------------------------------------|
| PCB 18                    |                              | 355±76 <sup>b</sup>                | 4.5±0.8°                              | 407±14 <sup>b</sup>                 |
| PCB 28                    | $62.8 \pm 1.2^{b}$           | 519±23 <sup>b</sup>                | $8.6 \pm 1.1^{c}$                     | $402 \pm 12^{b}$                    |
| PCB 44                    |                              | 417±31 <sup>b</sup>                |                                       | $408 \pm 21^{b}$                    |
| PCB 49                    |                              | $426\pm8^{\mathrm{b}}$             |                                       | 414±15 <sup>b</sup>                 |
| PCB 52                    | $7.65 \pm 0.23^{c}$          | $424 \pm 46^{b}$                   |                                       | $401 \pm 14^{b}$                    |
| PCB 66                    | $32.8 \pm 1.1^{b}$           | $429 \pm 39^{b}$                   | $5.6 \pm 1.7^{c}$                     | <b>414</b> ± <b>6</b> <sup>b</sup>  |
| PCB 74                    | $149\pm2^{\mathrm{b}}$       | 562±8 <sup>b</sup>                 | $13.8 \pm 0.1^{b}$                    | <b>414</b> ± <b>46</b> <sup>b</sup> |
| PCB 77                    | $0.673\pm0.047^{c}$          | $2.71\pm0.14^{c}$                  |                                       |                                     |
| PCB 81                    | $0.092\pm0.006^{c}$          | $0.630 \pm 0.028^{c}$              |                                       |                                     |
| PCB 87                    |                              | 455±52 <sup>b</sup>                |                                       | $396\pm23^{b}$                      |
| PCB 99                    | 139±5 <sup>b</sup>           | 552±8 <sup>b</sup>                 | $11.6 \pm 0.6^{c}$                    | $385 \pm 28^{b}$                    |
| PCB 101                   |                              | $441 \pm 29^{b}$                   |                                       | $409 \pm 27^{b}$                    |
| PCB 105                   | $47.5 \pm 6.7^{b}$           | $478\pm6^{\mathrm{b}}$             |                                       | 419±31 <sup>b</sup>                 |
| PCB 110                   | $12.3 \pm 0.4^{b}$           | 459±36 <sup>b</sup>                |                                       | $397 \pm 26^{b}$                    |
| PCB 114                   |                              | $90.5 \pm 7.4^{c}$                 |                                       | $46.6 \pm 8.4^{b}$                  |
| PCB 118                   | 213±4 <sup>b</sup>           | 594±34 <sup>b</sup>                | 18.9±1.2 <sup>b</sup>                 | 412±35 <sup>b</sup>                 |
| PCB 123                   |                              | 67.9±3.1°                          |                                       | $52.5 \pm 4.8^{b}$                  |
| PCB 126                   | $0.79\pm0.15^{b}$            | $10.4\pm2.5^{b}$                   |                                       | 7. 8±1.1 <sup>b</sup>               |
| PCB 128                   |                              | $400 \pm 29^{\rm b}$               |                                       | $420 \pm 17^{b}$                    |
| PCB 138                   | $319 \pm 27^{b}$             | $646 \pm 17^{\rm b}$               | $36.9 \pm 9.0^{b}$                    | 473±54 <sup>b</sup>                 |
| PCB 146                   | $55.6 \pm 0.8^{b}$           | 495±25 <sup>b</sup>                | $7.17\pm0.26^{c}$                     | $378 \pm 26^{b}$                    |
| PCB 149                   | $7.40\pm0.60^{c}$            | 409±3 <sup>b</sup>                 |                                       | $373 \pm 22^{b}$                    |
| PCB 151                   |                              | 418±11 <sup>b</sup>                |                                       | 381±18 <sup>b</sup>                 |
| PCB 153                   | 478±6 <sup>b</sup>           | 980±6 <sup>b</sup>                 | $58.2 \pm 0.9^{b}$                    | 457±36 <sup>b</sup>                 |
| PCB 156                   | $60.7 \pm 4.0^{\text{b}}$    | 505±7 <sup>b</sup>                 | $8.24 \pm 0.57^{b}$                   | 418±19 <sup>b</sup>                 |
| PCB 157                   | $14.7 \pm 0.6^{b}$           | $467 \pm 37^{\rm b}$               |                                       | 420±42 <sup>b</sup>                 |
| PCB 158                   |                              | <b>408</b> ± <b>8</b> <sup>b</sup> |                                       | 365±48 <sup>b</sup>                 |
| PCB 167                   | 15.9±2.2 <sup>b</sup>        | 467±59 <sup>b</sup>                |                                       | 403±26 <sup>b</sup>                 |
| PCB 169                   | $0.216\pm0.020^{b}$          | $9.3\pm2.0^{b}$                    |                                       | 8.10±0.33 <sup>b</sup>              |
| PCB 170                   | $100 \pm 18^{b}$             | 505±69 <sup>b</sup>                | $16.2 \pm 2.0^{b}$                    | 422±23 <sup>b</sup>                 |
| PCB 172                   | <b>14.8±0.6</b> <sup>b</sup> | 443±34 <sup>b</sup>                |                                       | 395±28 <sup>b</sup>                 |
| PCB 177                   | 23.7±1.2 <sup>b</sup>        | <b>450</b> ±6 <sup>b</sup>         |                                       | 392±9 <sup>b</sup>                  |
| PCB 178                   | 22.2±0.5 <sup>b</sup>        | 449±15 <sup>b</sup>                | $3.61\pm0.32^{c}$                     | 386±15 <sup>b</sup>                 |
| PCB 180                   | 234±8 <sup>b</sup>           | 678±58 <sup>b</sup>                |                                       | 459±49 <sup>b</sup>                 |
| PCB 180/193               |                              |                                    | <b>54.5</b> ± <b>0.5</b> <sup>b</sup> |                                     |
| PCB 183                   | 37.1±5.2 <sup>b</sup>        | 446±22 <sup>b</sup>                | $5.77 \pm 0.35^{\circ}$               | <b>407</b> ± <b>36</b> <sup>b</sup> |
| PCB 187                   | 96±17 <sup>b</sup>           | 517±28 <sup>b</sup>                | 15.5±0.5 <sup>b</sup>                 | 411±38 <sup>b</sup>                 |
| PCB 189                   |                              | 432±39 <sup>b</sup>                |                                       | 402±31 <sup>b</sup>                 |
| PCB 194                   | $50 \pm 16^{b}$              | 484±6 <sup>b</sup>                 | 11.9±0.3 <sup>b</sup>                 | 387±20 <sup>b</sup>                 |
| PCB 195                   | 11.9±1.6°                    | 459±19 <sup>b</sup>                |                                       | 385±25 <sup>b</sup>                 |
| PCB 196                   | 11.5 – 110                   | 107 = 17                           |                                       | 397±2 <sup>b</sup>                  |
| PCB 196/203               |                              | 870±31 <sup>b</sup>                | $10.9 \pm 0.7^{c}$                    | 0,7,-2                              |
| PCB 199                   |                              | 469±17 <sup>b</sup>                | 11.2±0.8°                             | 362±15 <sup>b</sup>                 |
| PCB 201                   |                              |                                    | 11.2-0.0                              | $397\pm7^{\rm b}$                   |
| PCB 203                   |                              |                                    |                                       | 398±49 <sup>b</sup>                 |
| PCB 206                   |                              | 463±12 <sup>b</sup>                | $6.98\pm0.37^{c}$                     | 366±8 <sup>b</sup>                  |
| PCB 209                   |                              | 441±26 <sup>b</sup>                | $3.33\pm0.59^{c}$                     | 338±16 <sup>b</sup>                 |
| 100 207                   |                              | 771-20                             | 5.55±0.57                             | 220-10                              |



Table 2 Certified, Reference, and *Information* mass fractions (ng/kg milk or reconstituted serum) of selected pesticides in SRM 1953 organic contaminants in non-fortified human milk, SRM 1954 organic contaminants in fortified human milk, SRM 1957 organic contaminants in non-fortified human serum, and SRM 1958 organic contaminants in fortified human serum

<sup>a</sup>Values are weighted means of the mass fractions determined by two to five analytical methods [29]. The uncertainty listed with each value is an expanded uncertainty about the mean [29, 35], with coverage factor, *k*=2, calculated by combining a pooled within method variance with a between method variance [29, 36] following the ISO Guide [37, 38]

bValues are the means of results using one analytical technique. The expanded uncertainty, U, is calculated as  $U=ku_c$ , where  $u_c$  is one standard deviation of the analyte mean, and the coverage factor, k, is determined from the Student's t-distribution corresponding to the associated degrees of freedom and 95 % confidence level for each analyte

| Pesticide                                                           | SRM 1953                             | SRM 1954                   | SRM 1957            | SRM 1958                     |
|---------------------------------------------------------------------|--------------------------------------|----------------------------|---------------------|------------------------------|
| 2,3,4,5-Tetrabromophenol                                            |                                      |                            | 120                 | 150                          |
| 2,3,4,6-Tetrachlorophenol                                           |                                      |                            | 99                  | 98                           |
| Pentachlorophenol                                                   |                                      |                            | $2710 \pm 880^a$    | $2780 \pm 550^{a}$           |
| Pentabromophenol                                                    |                                      |                            |                     | $445\!\pm\!93^a$             |
| Pentachloronitrobenzene                                             |                                      |                            |                     | 480                          |
| Pentachlorobenzene                                                  |                                      | $362{\pm}24^b$             |                     |                              |
| Hexachlorobenzene                                                   | $261\pm30^{a}$                       | $671 \pm 48^{a}$           | $29.7 \pm 3.5^{a}$  | $442 \pm 46^{a}$             |
| Octachlorostyrene                                                   |                                      | $374{\pm}28^b$             |                     | $359\!\pm\!19^a$             |
| α-НСН                                                               |                                      |                            |                     | $260\!\pm\!44^a$             |
| β-НСН                                                               | 610±4 <sup>a</sup>                   | <b>820</b> ±6 <sup>a</sup> | $31.3 \pm 6.0^a$    | $278\!\pm\!47^a$             |
| ү-НСН                                                               |                                      | 587±54 <sup>a</sup>        |                     | $315\!\pm\!43^a$             |
| Oxychlordane                                                        | $612\!\pm\!68^{\mathrm{a}}$          | $1060 \pm 150^{a}$         |                     | $226\!\pm\!43^a$             |
| cis-Chlordane                                                       |                                      | 368±9 <sup>a</sup>         |                     | $412\pm25^{a}$               |
| trans-Chlordane                                                     |                                      | $378{\pm}7^{\mathrm{a}}$   |                     | 412±3 <sup>a</sup>           |
| cis-Nonachlor                                                       | $126\!\pm\!6^{\mathrm{a}}$           | 499±5 <sup>a</sup>         |                     | 425±23 <sup>a</sup>          |
| trans-Nonachlor                                                     | $1240\!\pm\!70^{\mathrm{a}}$         | $1690 \pm 110^{a}$         | $58.3 \pm 1.9^{a}$  | $469 \!\pm\! 6^{\mathrm{a}}$ |
| Mirex                                                               | $68 \pm 25^{a}$                      | 512±7 <sup>a</sup>         |                     | $384\!\pm\!59^a$             |
| 2,4'-DDE                                                            |                                      | 403±6 <sup>a</sup>         |                     | $450\!\pm\!20^a$             |
| 4,4'-DDE                                                            | $7430 \!\pm\! 420^{\mathrm{a}}$      | $8150\pm200^{a}$           | 921±76 <sup>a</sup> | $1250\!\pm\!130^{\rm a}$     |
| 2,4'-DDD                                                            |                                      | $421 \pm 7^{a}$            |                     | $347\pm46^{a}$               |
| 4,4'-DDD                                                            |                                      | <b>422</b> ±8 <sup>a</sup> |                     | $416\!\pm\!17^{\mathrm{a}}$  |
| 2,4'-DDT                                                            | $18.2 \pm 1.2^{b}$                   | 443±9 <sup>a</sup>         |                     | 313±43 <sup>a</sup>          |
| 4,4'-DDT                                                            | $\textbf{229}\!\pm\!\textbf{16}^{a}$ | 691±37 <sup>a</sup>        |                     | $293\!\pm\!12^{\mathrm{a}}$  |
| 2-endo,3-exo,5-endo,6-exo,8,8,10,10-octachlorobornane (Parlar 26)   |                                      |                            | 340                 |                              |
| 2-endo,3-exo,5-endo,6-exo,8,8,9,10,10-nonachlorobornane (Parlar 50) |                                      |                            | 230                 |                              |
| 2,2,5,5,8,9,9,10,10-nonachlorobornane (Parlar 62)                   |                                      |                            | 320                 |                              |

non-fortified serum (SRM 1957) has equivalent or slightly higher mass fractions of the higher chlorinated PCB congeners, particularly PCB 180 (see Electronic Supplementary Material Figure S2). PCB 153/132 is the predominant PCB congener in both samples contributing almost 30 % of the total PCB composition (see Figure S2). PCB 153 and PCB 132 are separable on a number of GC columns; however, since a number of the methods for the certification of these SRMS relied on a 5 % phenyl methylpolysiloxane phase which does not separate PCB 132 from PCB 153, the certified value is for the combined PCB 153/132. On a lipid basis the non-fortified milk sample has higher mass fractions for most of the chlorinated pesticides and some of the PBDE congeners, particularly PBDE 100 and PBDE 153. For the PBDE congeners, PBDE 47 is the most abundant in both the milk and serum samples, contributing approximately 50 % of the total PBDE content (see Electronic Supplementary Material Figure S3). In comparing these samples, however, note that the milk and serum pools were not collected from the same individuals.

The goal in preparing the fortified serum SRM was to provide a material in which compounds not detectable in the average natural level found in the US would be present at detectable levels. These compounds may be present at higher concentrations in a sub-segment of the population with higher exposure; therefore reference materials with detectable levels are needed. Based on the certified values for the two serum SRMs, the PCB, pesticide, and PBDE mass fractions in the fortified material range from approximately 10 to 30 times the mass fractions in the non-fortified (natural) material with the exception of the highest level contaminants such as 4,4'-DDE and PBDE 47 which are only 1 to 3 times greater in SRM 1958. The two milk SRMs differ by only a factor of 2 to 4 in mass fraction levels for many of the pesticides and PBDEs (the predominant compounds), whereas the PCBs generally differ by a factor of 10 to 20. The fortified materials, SRM 1954 and SRM 1958, were both spiked with the same spiking solution (see Table S1). When the measured mass fraction differences between the fortified and non-fortified milk samples are compared to the difference between the fortified and non-fortified serum samples, the differences are similar and consistent. Based on the homogeneity analyses, the spiking solution was well mixed within the samples prior to bottling.

In addition to the compounds characterized in this study, SRM 1954 and SRM 1958 were fortified with compounds



**Table 3 Certified** and Reference mass fractions (ng/kg milk or reconstituted serum) of selected PBDE congeners and PBB 153 in SRM 1953 organic contaminants in non-fortified human milk, SRM 1954 organic contaminants in fortified human milk, SRM 1957 organic contaminants in non-fortified human serum, and SRM 1958 organic contaminants in fortified human serum

| PBDE congener <sup>a</sup> | SRM 1953                    | SRM 1954                      | SRM 1957              | SRM 1958                  |
|----------------------------|-----------------------------|-------------------------------|-----------------------|---------------------------|
| PBDE 17                    | 16.8±0.1 <sup>b</sup>       | 401±50 <sup>b</sup>           | 4.2±1.5 <sup>b</sup>  | 458±32 <sup>b</sup>       |
| PBDE 28                    | $137\!\pm\!18^b$            | $576 \pm 7^{b}$               | $20.0\!\pm\!2.4^b$    | $462 \pm 19^{b}$          |
| PBDE 47                    | $2240 \pm 30^{b}$           | $2590\!\pm\!80^{\mathrm{b}}$  | $268 \pm 14^{b}$      | $651 \pm 29^{b}$          |
| PBDE 66                    | $16.2 \pm 3.4^{b}$          | $414{\pm}10^{\mathrm{b}}$     | $6.70\!\pm\!0.13^{b}$ | $440{\pm}41^{\rm b}$      |
| PBDE 85                    | $49.2 \pm 8.1^{b}$          | $475{\pm}22^{\mathrm{b}}$     | $8.2 \pm 1.9^{b}$     | $475 \pm 39^{b}$          |
| PBDE 99                    | $341{\pm}8^{\rm b}$         | $738{\pm}47^{\mathrm{b}}$     | $76.0 \pm 3.8^{b}$    | $492 \pm 15^{\mathrm{b}}$ |
| PBDE 100                   | $822\!\pm\!15^{\mathrm{b}}$ | $1280\!\pm\!40^{\mathrm{b}}$  | $49.7 \pm 2.7^{b}$    | $475{\pm}27^{\mathrm{b}}$ |
| PBDE 153                   | $971 \pm 23^{\mathrm{b}}$   | $1440\!\pm\!100^{\mathrm{b}}$ | $61.0 \pm 3.2^{b}$    | $455{\pm}54^{\rm b}$      |
| PBDE 154                   | $42.7 \pm 2.3^{b}$          | $470\!\pm\!13^{\rm b}$        | $7.0\!\pm\!1.0^{c}$   | $441 \pm 39^{b}$          |
| PBDE 183                   |                             | $507 \pm 43^{b}$              |                       | $453 \pm 42^{b}$          |
| PBDE 203                   |                             | $494\!\pm\!44^c$              |                       |                           |
| PBDE 206                   |                             | $820\!\pm\!280^b$             |                       | $426\pm4^{\mathrm{b}}$    |
| PBDE 209                   |                             | $423\!\pm\!24^{\rm b}$        |                       | $417\pm5^{b}$             |
| PBB 153                    | $36.5 \pm 0.6^{b}$          | 476±9 <sup>b</sup>            | $15.5 \pm 0.1^{b}$    | $421\!\pm\!13^{\rm b}$    |

<sup>&</sup>lt;sup>a</sup> PBDE congeners and PBB 153 are numbered according to IUPAC rules

anticipated to be of interest in future monitoring studies including additional brominated flame retardants, polybrominated dibenzodioxins, polybrominated dibenzofurans, mixed chlorinated and brominated dibenzodioxins and dibenzofurans, polychlorinated naphthalenes, and hydroxylated PCBs (see Table S1). As new flame retardants are introduced into products, they will need to be monitored in human samples, as well as the combustion products from the incineration of the products, which produce the brominated and mixed chlorinated/brominated dioxins and furans. Hydroxylated metabolites of the contaminants, such as the hydroxylated PCBs are toxicologically more active than their parent compounds, although found at lower levels in the serum and milk. With these additional compounds added to the serum and milk pools, these SRMs will be valuable for future method development. The serum materials have also been characterized for selected perfluorinated compounds (PFCs) as described in Keller et al. [30]. (PFCs were

Table 4 Reference mass fractions (pg/kg milk or reconstituted serum) of selected PCDD and PCDF congeners in SRM 1953 organic contaminants in non-fortified human milk, SRM 1954 organic contaminants in fortified human milk, SRM 1957 organic contaminants in non-fortified human serum, and SRM 1958 organic contaminants in fortified human serum

| Dioxin or furan congener <sup>a</sup> | SRM 1953              | SRM 1954             | SRM 1957             | SRM 1958              |
|---------------------------------------|-----------------------|----------------------|----------------------|-----------------------|
| 2,3,7,8-TCDD                          | 55.2±4.6 <sup>a</sup> | 162±32 <sup>a</sup>  |                      | 94.2±7.8 <sup>a</sup> |
| 1,2,3,7,8-PCDD                        | $123\!\pm\!11^a$      | $240\!\pm\!27^a$     | $17.3\!\pm\!2.9^{a}$ | $106\!\pm\!12^a$      |
| 1,2,3,4,7,8-HxCDD                     | $77.7 \pm 9.0^{a}$    | $180{\pm}21^a$       | $11.8\!\pm\!1.7^a$   | $95.3\!\pm\!9.9^a$    |
| 1,2,3,6,7,8-HxCDD                     | $501\!\pm\!91^a$      | $890\!\pm\!240^a$    | $82\!\pm\!18^a$      | $340\!\pm\!60^a$      |
| 1,2,3,7,8,9-HxCDD                     | $100\!\pm\!11^a$      | $206{\pm}33^a$       | $18.7{\pm}0.3^a$     | $99.6{\pm}7.2^a$      |
| 1,2,3,4,6,7,8-HpCDD                   | $480\!\pm\!160^{a}$   | $1080\!\pm\!410^a$   | $104{\pm}36^a$       | $565{\pm}90^a$        |
| OCDD                                  | $2240{\pm}360^a$      | $4890\!\pm\!1430^a$  | $716{\pm}81^a$       | $2570\!\pm\!280^a$    |
| 2,3,7,8-TCDF                          | $16.4{\pm}0.6^a$      | $125\!\pm\!13^a$     |                      | $104{\pm}3^a$         |
| 1,2,3,7,8-PCDF                        |                       | $132{\pm}29^a$       |                      | $98\!\pm\!15^a$       |
| 2,3,4,7,8-PCDF                        | $127\!\pm\!11^a$      | $344{\pm}33^a$       | $16.4{\pm}2.5^a$     | $199{\pm}30^a$        |
| 1,2,3,4,7,8-HxCDF                     | $72\!\pm\!10^a$       | $171\!\pm\!22^a$     | $17.0{\pm}3.1^a$     | $95.9{\pm}9.3^a$      |
| 1,2,3,6,7,8- HxCDF                    | $66.0 \pm 1.3^a$      | $183\!\pm\!20^a$     | $14.1{\pm}2.5^a$     | $102\!\pm\!11^a$      |
| 1,2,3,7,8,9- HxCDF                    | $31.6{\pm}0.8^a$      | $126{\pm}9^a$        |                      | $94\pm11^a$           |
| 2,3,4,6,7,8- HxCDF                    |                       | $1080\!\pm\!160^{a}$ |                      | $900 \pm 110^{a}$     |
| 1,2,3,4,6,7,8-HpCDF                   | $100\!\pm\!12^a$      | $405{\pm}66^a$       | $39.8\!\pm\!6.2^a$   | $289{\pm}9^a$         |
| 1,2,3,4,7,8,9-HpCDF                   |                       |                      |                      | $84\!\pm\!12^a$       |
| OCDF                                  | $11.6 \pm 5.0^{b}$    | $94.3 \pm 12.7^{b}$  |                      | $83.3\!\pm\!1.7^{a}$  |

<sup>&</sup>lt;sup>a</sup> Values are weighted means of the mass fractions determined by two to three analytical methods [29]. The uncertainty listed with each value is an expanded uncertainty about the mean [29, 35], with coverage factor, k=2, calculated by combining a pooled within method variance with a between method variance [29, 36] following the ISO Guide [37, 38]

not included in the spiking solution.) Nutrient elements, including calcium, copper, iron, magnesium, manganese, phosphorous, potassium, sodium, and zinc, and the toxic metal mercury, were determined in the milk samples, and the results are provided on the Certificates on Analysis as reference mass fractions [18, 20, 31].

SRM 1957 was developed in part to address the need for a serum material representing more contemporary contaminant levels in the U.S. population than the previous material, SRM 1589a. Using biomonitoring data, Sjödin et al. [32] demonstrated the decrease in legacy chemicals (no longer in use in the U.S.) such as PCBs and the increase in more recently used chemicals such as PBDEs in serum collected from 1984 to 2002. Comparison of the mass fractions of contaminants in SRM 1589a (collected in 1996) and SRM 1957 (collected in 2004) indicates similar trends as shown in Fig. 1. For example, PCB congener values in SRM 1957 are




<sup>&</sup>lt;sup>b</sup> Values are weighted means of the mass fractions determined by two to five analytical methods [29]. The uncertainty listed with each value is an expanded uncertainty about the mean [29, 35], with coverage factor, *k*=2 (except *k*=4 for PBDE 206 in SRM 1958), calculated by combining a pooled within method variance with a between method variance [29, 36] following the ISO Guide [37, 38]

<sup>&</sup>lt;sup>c</sup> Values are the means of results using one analytical technique. The expanded uncertainty, U, is calculated as  $U=ku_c$ , where  $u_c$  is one standard deviation of the analyte mean, and the coverage factor, k, is determined from the Student's t-distribution corresponding to the associated degrees of freedom and 95 % confidence level for each analyte

<sup>&</sup>lt;sup>b</sup> Values are the means of results using one analytical technique. The expanded uncertainty, U, is calculated as  $U=ku_c$ , where  $u_c$  is one standard deviation of the analyte mean, and the coverage factor, k, is determined from the Student's t-distribution corresponding to the associated degrees of freedom and 95 % confidence level for each analyte

Fig. 1 Comparison of mass fractions for legacy chemicals (PCBs and 4,4'-DDE) and chemicals used more recently (PBDEs and PFCs) in serum collected in 1996 (SRM 1589a) and serum collected in 2004 (SRM 1957). Error bars represent the expanded uncertainties

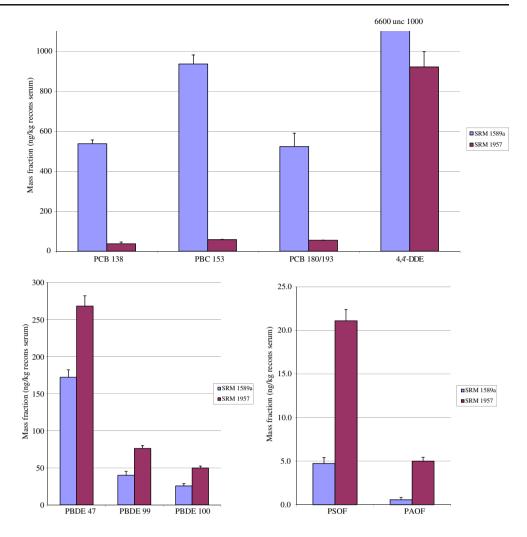

generally 10 to 15 times lower than in SRM 1589a. Similarly the mass fraction of 4,4'-DDE was 12 times lower in SRM 1957 compared to SRM 1589a. For the PCDD/PCDF congeners (not shown in Fig. 1), the mass fractions were typically a factor of four lower in SRM 1957 compared to SRM 1589a, e.g., 1,2,3,6,7,8-hexaCDD  $(416 \text{ pg/kg}\pm 12 \text{ pg/kg vs. } 84.7 \text{ pg/kg}\pm 8.1 \text{ pg/kg})$ . For the flame retardants such as PBDEs, the values in SRM 1957 have increased typically by a factor of two compared to SRM 1589a. Likewise for PFCs, another group of compounds used recently, mass fractions were a factor of 5 and 10 higher in SRM 1957 compared to SRM 1589a for the predominant PFCs [30]. In making the comparison of SRM 1589a to SRM 1957, however, it should be noted that SRM 1589a was prepared from serum collected from donor units that had been prescreened at NIST for the PCB content, and the units with high PCB content were selected (all coming from the Chicago area, preferably from donors who ate fish caught in the Great Lakes), whereas SRM 1957 was prepared from random donor units (i.e., no prescreening) from multiple collection sites.

#### Conclusion

These new human serum and milk SRMs represent human samples commonly used for the monitoring of environmental contaminants. SRMs 1953 and 1954 are the first human milk certified reference materials available worldwide that have been characterized for organic contaminants. It should be noted that although there are some environmental contaminants present in human milk, the nutritional benefits of breastfeeding to the infant likely outweigh any harmful effects (for example, see http://www.cdc.gov/ breastfeeding/). SRM 1957 represents the most extensively characterized natural human serum certified reference material for persistent organic contaminants available from any source worldwide. The fortified milk and serum SRMs, with certified and reference values for 91 and 92 PCBs, pesticides, PBDEs, and dioxins/furans, respectively, are valuable reference materials for contaminants currently monitored, as well as for contaminant classes that may be monitored extensively in the future. These materials will find extensive use in the biomonitoring community and will aid in the



comparison of data among biomonitoring studies performed in different years and by different laboratories.

Acknowledgments We would like to thank Kim Updegrove (Mothers' Milk Bank at Austin, Austin, TX) and Pauline Sakomoto (Mothers' Milk Bank, San Jose, CA) for assistance with acquisition of the milk. Laboratories participating in an interlaboratory study included Institut National de Santé Publique du Québec (Québec, Canada) coordinated by Pierre Dumas and Stockholm University (Stockholm, Sweden) coordinated by Maria Athanasiadou.

**Disclaimer** Certain commercial equipment, instruments, or materials are identified in this paper to specify adequately the experimental procedure. Such identification does not imply recommendation or endorsement by the National Institute of Standards and Technology, nor does it imply that the materials or equipment identified are necessarily the best available for the purpose.

## References

- Sjödin A, Wong L-Y, Johes RS, Park A, Ahang Y, Hodge C, Dipietro E, McClure C, Turner W, Needham LL, Patterson DG Jr (2008) Environ Sci Technol 42:1377–1384
- Patterson DG Jr, Wong L-Y, Turner WE, Caudill SP, Dipietro ES, McClure PC, Cash TP, Osterloh JD, Pirkle JL, Sampson EJ, Needham LL (2009) Environ Sci Technol 43:1211–1218
- Sandanger TM, Sinotte M, Dumas P, Marchand M, Sandau CD, Pereg D, Bérubé S, Brisson J, Ayotte P (2007) Environ Health Perspect 115:1429–1434
- Roosens L, Abdallah MA-E, Harrad S, Neels H, Covaci A (2009) Environ Sci Technol 43:3535–3541
- Inoue K, Harada K, Takenaka K, Uehara S, Kono M, Shimizu T, Takasuga T, Senthilkumar K, Yamashita F, Koizumi A (2006) Environ Health Perspect 114:1179–1185
- Kim B-H, Ikonomou MG, Lee S-J, Kim H-S, Chang Y-S (2005) Sci Total Environ 336:45–56
- Hagmar L, Wallin E, Vessby B, Jönsson BAG, Bergman A, Rylander L (2006) Chemosphere 64:1507–1513
- Gladen BC, Monaghan SC, Lukyanova EM, Hulchiy OP, Shkyryak-Nyzhnyk ZA, Sericano JL, Little RE (1999) Environ Health Perspect 107:459–462
- Paustenbach D, Galbraith D (2006) Environ Health Perspect 114:1143–1149
- Albertini R, Bird M, Doerrer N, Needham L, Robison S, Sheldon L, Zenick H (2006) Environ Health Perspect 114:1755–1762
- Patterson DG Jr, Turner WE, Caudill SP, Needham LL (2008) Chemosphere 73:5261–5277
- Certificate of Analysis, SRM 1589 Polychlorinated Biphenyls (as Aroclor 1260) in Human Serum (1985) NIST, Gaithersburg. https://www-s.nist.gov/srmors/certificates/view\_certGIF.cfm? certificate=1589. Accessed 7 Sept 2012
- Certificate of Analysis, SRM 1589a PCBs, Pesticides, and Dioxins/Furans in Human Serum (2000) NIST, Gaithersburg. https://www-s.nist.gov/srmors/certificates/archive/1589a% 20August%209,%202000.pdf. Accessed 7 Sept 2012
- Certificate of Analysis, SRM 1589a PCBs, Pesticides, PBDEs, and Dioxins/Furans in Human Serum. (2006) NIST, Gaithersburg. https://www-s.nist.gov/srmors/certificates/view\_certGIF.cfm? certificate=1589a. Accessed 7 Sept 2012
- Schantz MM, Keller JM, Leigh S, Patterson DG, Sharpless KE, Sjodin A, Stapleton HM, Swarthout R, Turner W, Wise SA (2007) Anal Bioanal Chem 389:1201–1208
- 16. Hites RA (2004) Environ Sci Technol 38:945-956

- Bernert JT, Turner WE, Patterson DG Jr, Needham LL (2007) Chemosphere 68:824–831
- Certificate of Analysis, SRM 1953 Organic Contaminants in Non-Fortified Human Milk (2009) NIST, Gaithersburg. https://www-s.nist. gov/srmors/certificates/view\_certGIF.cfm?certificate=1953. Accessed 7 Sept 2012
- Certificate of Analysis, SRM 1957 Organic Contaminants in Non-Fortified Human Serum (2009) NIST, Gaithersburg. https://www-s. nist.gov/srmors/view\_cert.cfm?srm=1957. Accessed 7 Sept 2012
- Certificate of Analysis, SRM 1954 Organic Contaminants in Fortified Human Milk (2009) NIST, Gaithersburg. https://www-s. nist.gov/srmors/certificates/view\_certGIF.cfm?certificate=1954. Accessed 7 Sept 2012
- Certificate of Analysis, SRM 1958 Organic Contaminants in Fortified Human Serum (2009) NIST, Gaithersburg. https://www-s.nist.gov/ srmors/certificates/view\_certGIF.cfm?certificate=1958. Accessed 7 Sept 2012
- Keller JM, Swarthout RF, Carlson BKR, Yordy J, Guichard A, Schantz MM, Kucklick JR (2009) Anal Bioanal Chem 393: 747, 760
- Patterson DG Jr, Turner WE (1997) CLIA document.
   Environmental Health Laboratory, Center for Environmental Health, Centers for Disease Control and Prevention, Atlanta
- Sjödin A, Jones RS, Lapeza CR, Focant J-F, McGahee EE III, Patterson DG Jr (2004) Anal Chem 76:1921–1927
- Sjödin A, McGahee EE III, Focant J-F, Jones RS, Lapeza CR, Zhang Y, Patterson DG Jr (2004) Anal Chem 76:4508–4514
- Focant J-F, Eppe G, Massart A-C, Scholl G, Pirard C, De Pauw E (2006) J Chromatogr A 1130:97–107
- Wise SA, Poster DL, Kucklick JR, Keller JM, Vander Pol SS, Sander LC, Schantz MM (2006) Anal Bioanal Chem 386: 1153–1190
- 28. May W, Parris R, Beck C, Fassett J, Greenberg R, Guenther F, Kramer G, Wise S, Gills T, Colbert J, Gettings R, MacDonald B (2000) NIST Special Publication 260–136, U.S. Government Printing Office: Gaithersburg. http://ts.nist.gov/MeasurementServices/ReferenceMaterials/PUBLICATIONS.cfm. Accessed 7 Sept 2012
- 29. Dersimonian R, Laird N (1986) Control Clin Trials 7:177-188
- Keller JM, Calafat AM, Kato K, Ellefson ME, Reagen WK, Strynar M, O'Connell S, Butt CM, Mabury S, Small J, Muir D, Leigh SD, Schantz MM (2010) Anal Bioanal Chem 307:439–451
- Sharpless KE, Lindstrom RM, Nelson BC, Phinney KW, Rimmer CA, Sander LC, Schantz MM, Spatz RO, Thomas JB, Turk GC, Wise SA, Wood LJ, Yen JH (2010) J AOAC Int 93:1262–1274
- Sjödin A, Jones RS, Focant JF, Lapeza C, Wang RY, McGahee EE, Zhang YL, Turner WE, Slazyk B, Needham LL, Patterson DG Jr (2004) Environ Health Perspec 112:654–658
- 33. Ballschmiter K, Zell M (1980) Fresenius Z Anal Chem 302:20-31
- 34. Schulte E, Malisch R (1983) Fresenius Z Anal Chem 314:545-551
- 35. Rukhin AL (2009) Metrologia 46:323-331
- 36. Horn RA, Horn SA, Duncan DB (1975) J Am Stat Assoc 70: 380–385
- 37. JCGM 100:2008; Joint Committee for Guides in Metrology (2008). available at http://www.bipm.org/utils/common/documents/jcgm/JCGM\_100\_2008\_E.pdf. Accessed 7 Sept 2012; see also Taylor, B.N.; Kuyatt, C.E.; Guidelines for Evaluating and Expressing the Uncertainty of NIST Measurement Results; (1994) NIST Technical Note 1297; U.S. Government Printing Office: Washington, DC. available at http://www.nist.gov/physlab/pubs/index.cfm. Accessed 7 Sept 2012
- JCGM 101:2008; Joint Committee for Guides in Metrology (BIPM, IEC, IFCC, ILAC, ISO, IUPAC, IUPAP and OIML) (2008) International Bureau of Weights and Measures (BIPM), Sèvres, France. available at http://www.bipm.org/utils/common/ documents/jcgm/JCGM\_101\_2008\_E.pdf. Accessed 7 Sept 2012

